
80K SLIMMING MACHINE


< Instruction Manual >



In order to use this product correctly, please read the manual carefully and keep it properly.

Part name

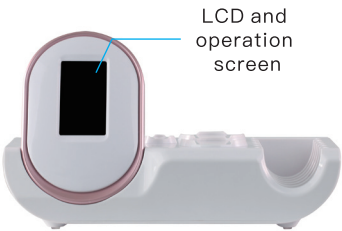

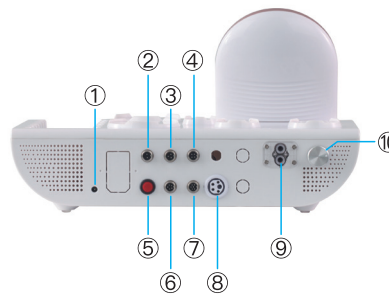

① Power interface

② 80K interface

③ 3–pole RF interface

④ 4–pole RF interface

⑤ Switch

⑥ Laser board interface

⑦ 6–pole RF interface

⑧ Negative pressure interface

⑨ Filter

⑩ Negative pressure regulation



3 RF Head



4 RF Head




6 RF Head



Negative pressure massage Head




80K Fat Burn Head

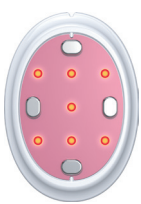


Laser &EMS

Usage method

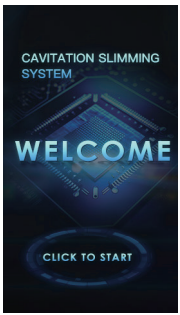





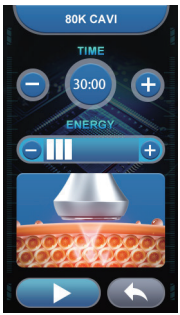

1. Connect the functional head to the corresponding socket and connect the power cord.
2. Press the power switch, the indicator light and display screen will light up. At this point, click on the touch screen to start.
3. Choose the job head you want.
4. After selection, adjust the required working time and working energy, click “ ▶ ” button to start working, click“ ⏸ ”button to suspend use, or click“ ↩ ”to return to the working head selection interface.
5. After use, please turn off the power.

Function of Different Head

3 RF Head:Suitable for 1–2 Energy



When using, select the 3–pole RF head and enter the function panel.
TIME: The default is 30 minutes, adjustable from 0 to 60 minutes. When it is 0, the instrument will stop working.
ENERGY: Adjust the RF energy of the working head. The larger the number, the greater the energy. The 10th gear can be adjusted.

Use for the face or eyes

When use on the skin, it can firms and lifts the skin, reduces winkle.Improvement eyes circulation, shrink pores, and puffy eyes.Eliminate wrinkles, dark circles and bags under the eyes.Improve skin texture and eliminate sagging eyes.Rejuvenates the skin and improves eye contour.

4 RF Head:Suitable for 2–3 Energy



When using, select the 4–pole RF head and enter the function panel.
TIME: The default is 30 minutes, adjustable from 0 to 60 minutes. When it is 0, the instrument will stop working.
ENERGY: Adjust the RF energy of the working head. The larger the number, the greater the energy. The 10th gear can be adjusted.

Use for the face and neck

Firm skin, enhance skin elasticity, facial wrinkle elimination, skin tightening and lifting,reduce face swelling, reduce redness, relieve eye fatigue, remove dark circles, And eye bags, improve skin quality.

6 RF Head:Suitable for 2–5 Energy




When using, select the 6–pole RF head and enter the function panel.
TIME: The default is 30 minutes, adjustable from 0 to 60 minutes. When it is 0, the instrument will stop working.
ENERGY: Adjust the RF energy of the working head. The larger the number, the greater the energy. The 10th gear can be adjusted.

Use for the back and waist of Body, thighs and calves

Can tighten the waist, firm and help body shaping.Soften the muscles of the arms and thighs, give the whole body massage and tighten the arms and thighs.

80K Fat Burn Head

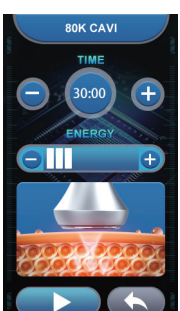

When using, select the 80K Fat Burn Head and enter the function panel.
TIME: The default is 30 minutes, adjustable from 0 to 60 minutes. When it is 0, the instrument will stop working.
ENERGY: Adjust the vibration energy of the working head. The larger the number, the greater the energy. The 10th gear can be adjusted.

Use for the body Shaping

It can dissolve fat, detoxify lymph, tighten skin, enhance skin elasticity, reach fat layer directly, have the function of directional RF output, accelerate the movement of fat cells to generate heat energy, and discharge the excess fat and toxin in the face out of the body through sweat glands, liver intestine circulation and lymph circulation, so as to dissolve fat.

Negative pressure massage Head



When using, select the Negative pressure massage Head and enter the function panel.
TIME: The default is 30 minutes, adjustable from 0 to 60 minutes. When it is 0, the instrument will stop working.
ENERGY: Adjust the RF energy of the working head. The larger the number, the greater the energy. The 10th gear can be adjusted.
SUCTION: Adjust the suction time of negative pressure. The larger the number is, the longer the suction time will be. It is invalid when it is 0.
RELEASE: Adjust the suction time of negative pressure. The larger the number is, the longer the suction time is. It is invalid when it is 0.
Note: The inspiratory and outgassing intervals are only available when the number is adjusted at the same time.

Body maintenance

Diversified level adjustment make your body to be care comfortably.Improve skin texture by promoting the regeneration of collagen and elasticity.Improve skin sagging caused by fat and skin aging.Help eliminate waste materials in the body and promote metabolism.Negative pressure massage function can promote fat cell volume reduction.Vacuum negative pressure massage can fully absorb various skin tissues between the multi–polar thermal melting radio frequency, so that the radio frequency energy can reach the deep skin tissue more effectively.At the same time, the

adsorption of negative pressure promotes blood circulation and metabolism. Let the muscle tissue be fully relaxed and relieve the fatigue of the skin and muscle tissue.

Laser & EMS



When using, select the Laser Lipolysis and enter the function panel.
TIME: The default is 30 minutes, adjustable from 0 to 60 minutes. When it is 0, the instrument will stop working.
Energy: Adjust the micro electric energy of the laser board, the larger the number, the greater the energy. 10 gears can be adjusted.
Light: Adjust the brightness of the laser board. The larger the number, the higher the brightness. 10 gears can be adjusted.


Use for the body Shaping

The Following information is provided for the correct utilization of IPL Skin treatment system.the Information includes not only the accident protection regulations the products comply with,but also the effective regarding proper use of the products.

Using Guidance

Thigh

(1) From bottom to top, push to groin to dredge lymph;
(2) Turning anticlockwise from bottom to top can help break down fat;
(3) Decompose greasy by pushing back and forth;
(4) Pull down from the knee to the groin to improve the curve.

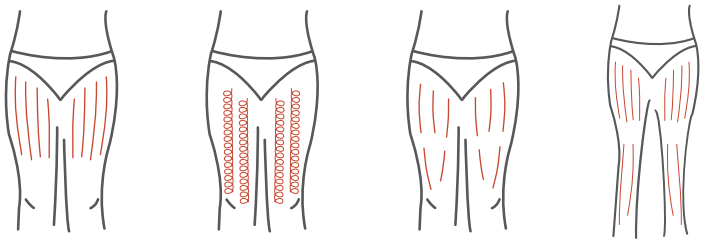


Hips

(1) Pull to the waist in the direction of lymph;
(2) Circle anticlockwise from bottom to top and move back and forth;
(3) According to the muscle texture, it can slide and compact the buttock skin, making the buttock plump and firm, presenting the perfect curve.



Back

(1) Push back and forth 2–3 times;
(2) Circle anticlockwise to stimulate the lymph nodes;
(3) On the back, turn anticlockwise and move up and down to help break down fat;
(4) Using lymphatic drainage technique to pull toxin to lymph nodes can help improve back curve.

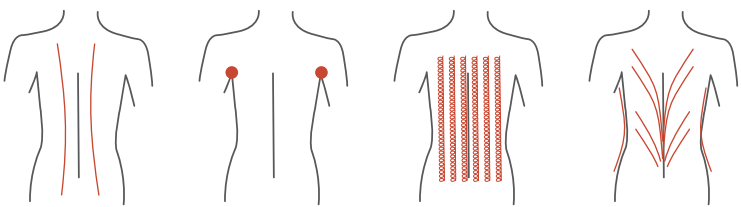


Arm

(1) From the inner elbow to the armpit, do lymphatic drainage;
(2) To stimulate lymph nodes in a circle;
(3) Break down excess fat by turning it anticlockwise orpushing it back and forth;
(4) Tighten skin from elbow to armpit.




Face


(1) Evenly apply the product on the face, use the face RF head, operate from the bottom to the top to lift the face, make the skin in a tense state, it is recommended to use for 5 minutes;
(2) Apply the product evenly to the face, use the RF head on the face, operate in the direction of the arrow, it is recommended to use it for 5–8 minutes.

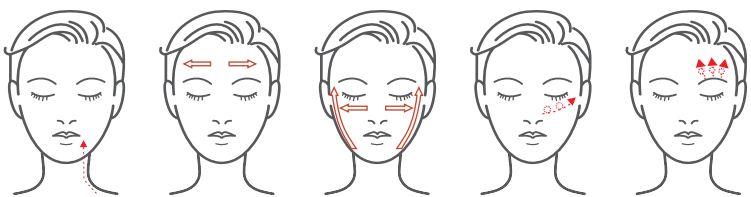

Abdomen

(1) Around the navel from small to large, clockwise circle to help large intestine peristalsis;
(2) From the small circle to the large circle, push it anticlockwise.

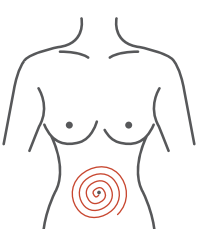

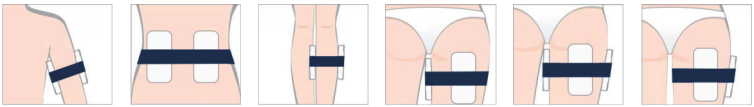

Upper arm | Lower abdomen/ back/sides | Calf | Thigh | Thigh | Thigh

Warm Tips for You

1. Do not use on your original and dry skin surface, it will burn your skin. Make sure you apply cream or oils to your skin then you can use.

2. When using, there must be an aqueous medium (essential oil/creams, etc.),otherwise it will damage to the machine Do not use on the same part for too long, and avoid body parts with bones.

3. Do not connect the power until probe is well installed.

4. The machine is not allowed to replace probe when in the charging state.

5. Please take off the metal jewelry when using the machine.

6. When the machine is used for 1 hour, please stop for 10 minutes to let the machine rest.

7. Make sure you turn off the power after using, and clean the probe keep good maintenance for next time use.

8. Please keep the probe clean every time, use a dry cloth to wipe it, do not use water to clean.

9. After using, clean you skin and apply hydrating mask or creams, oils after 10 minutes.

10. When using 80KHead, some people may experience tinnitus. This is normal because this sound wave exceeds the range of human sound waves.

Warranty & After Service

We promise that all machines undergo our rigorous inspection to ensure the best quality. For usage issues, we will use our professional skills to address your concerns.

If you have any questions or other questions about our products, please feel free to contact us at any time. We will always support you.

The host comes with a one-year free warranty and a three-year paid warranty. Accessories, power supplies, etc. are not covered by the warranty!

Made in China